# Operation Microscope OMS-90




# An innovative concept providing optimal observation and versatile operation



# Coaxial Illumination

The newly designed illumination system of the OMS-90 assures the finest optical performance. The optical system has been fitted with a new and innovative design that allows the surgeon to view a sharper image, with more clarity than ever before.

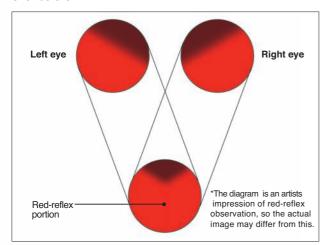

The illumination system of the OMS-90 produces a superior red-reflex image forassistance in surgical procedures. As shown in the diagram, theinnovative system design merges two improved field of views to produce an enhanced red-reflex andunparalleled image quality.

# Stain-proof Coated Objective Lens (F=200mm)

The newly developed stain-proof coating is applied to help retain a clear viewing area. The focal length (F=200mm) provides the standardized comfortable working distance.

# Parallel Binocular Tube

The parallel binocular tube provides a natural stereoscopic view.

# PD Adjustment Knob

The built-in PD adjustment feature offers easy aligment.

# Built-in Yellow Filter

The flip in/out yellow filter is built into the system to allow full protection against retina phototoxicity.






# Compact Base

The new sturdy and compact \$\phi740mm\$ (30") base is designed for optimal use of limited operating room floor space.

**OMS-90** 

Previous model

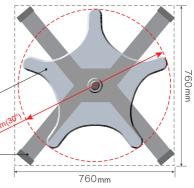

# Optional Accessories

# X-Y Translator

The X-Y translator adapter provides precise movement around the surgical field. A one-touch re-centering button is provided.



# Assistant's Microscope

The parallel binocular assistant's microscope provides easy positioning, a wide field of view and 30 degrees of rotation.



# Beam Splitter (80/20 or 50/50)

An 80/20 or a 50/50 beam splitter is available for the attachment of documentation devices.

# Oblique Illuminator

The oblique illuminator offers a supplemental light source for indirect illumination. The slit will rotate 30 degrees and adjust from 0 to 20mm in width.




# Accessory Adapter Model 2

The accessory adapter is required to attach the assistant's microscope and/or the oblique illuminator to the OMS-90.

# 1/2 TV Attachment

The 1/2 TV attachment incorporates a built-in beam splitter and lets the user connect a 1/2 inch C-mount camera to the OMS-90.

# **System Chart**

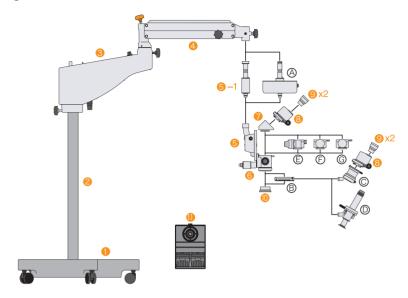

# **Specifications**

# Main microscope

Tidiii iiiidi dadapa	
Binocular tube	Parallel 45° inclined binocular tube
Objective lens	F=200mm
Eyepieces	12.5x
Magnification change type	Drum type with 5 changes
Total magnification	3.4x, 5.3x, 8.5x, 13.6x, 21.2x
Field of view	φ 66.3, φ 42.3, φ 26.3, φ 17.1, φ 1 0.3mm
X-Y Translator (option)	Range:±25mm

# Arm(counter-balanced)

Length	1st arm:400mm; 2nd arm:560mm
Rotational angle	1st arm:360°; 2nd arm:300°
Column base	
Base dimension	φ 740mm (30")
Power unit	
Primary voltage	AC100~120V or AC220~240V
Power frequency	50/60Hz
Maximum power consumption	280VA
Weight	
Weight	97kg
<u> </u>	<u> </u>

\*Subject to change in design and/or specifications without advanced notice. In compliance with the terms of the Export Administration Regulation of the United States of America, this product may not be available in some regions or countries.



IMPORTANT In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.





2 ----- Column

-----1st Arm

-----2nd Arm

----- Focusing Unit

6 −1····Joint Arm

----- Microscope Body

-----45° Inclined Adapter

8----- Parallel Binocular Tube

9-----12.5X Eyepiece

10------F=200mm Objective Lens

10-----Foot Switch

# < Optional Accessories >

(A)----- X-Y Translator

Accessory Adapter Model 2

O----- Assistant's Microscope

🔘----- Oblique Illuminator

(E)----- 1/2 TV Attachment

(F)---- Beam Splitter(80/20)

(G)----- Beam Splitter(50/50)







TOPCON CORPORATION

TOPCON CORPORATION



# **TOPCON** TOPCON CORPORATION

75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, Japan. Phone:3-3558-2523/2522 Fax:3-3960-4214 www.topcon.co.jp

TOPCON MEDICAL SYSTEMS, INC. 111 Bauer Drive, Oakland, NJ 07436, U.S.A. Phone:+1-201-599-5100 Fax:+1-201-599-5250 www.topconmedical.com

TOPCON CANADA INC.
110 Provencher Avenue, Boisbriand, QC J7G 1N1 CANADA
Phone:+1-450-430-7771 Fax:+1-450-430-6457 www.topco

# TOPCON EUROPE MEDICAL B.V.

Essebaan 11; 2908 LJ Capelle a/d IJssel; P.O.Boxt 2900 AC Capelle a/d IJssel; THE NETHERLANDS Phone:41 (0)10-4585045 E-mail: medical@topcon.nl; www.topcon.eu

Viale dell' Industria 60; 20037 Paderno Dugnano; (Milano), ITALY
Phone:+39-02-9186671 Fax:+39-02-91081091 E-mail: topconitaly@tiscali.it; www.topcon.it

# DANMARK OFFICE Praestemarksvej 25; 4000 Roskilde, DANMARK Phone: 45-46-327500 Fax: -45-46-327555 E-mail: topcon@topcondanmark.dk www.topcondanmark.dk

IRELAND OFFICE

Unit 276, Blanchardstown; Corporate Park 2 Ballycoolin Dublin 15, IRELAND Phone:+353-18975900 Fax:+353-18293915 E-mail: medical@topcon.ie; www.topcon.ie

# TOPCON S.A.R.L.

DANMARK OFFICE

rue de Paris 92585 Clichy, FRANCE ne:+33-(0)1-41069494 Fax:+33-(0)1-47390251 E-mail:topcon@topcon.fr; www.topcon.fr

TOPCON DEUTSCHLAND G.m.b.H.
Hanns-Marlin-Schleyer Strasse 41; D-47877 Willich, GERMANY
Phone: 49-(0)2154-8850 Fax: +49-(0)2154-885177 E-mail:med@topcon.de; www.topcon.de

TOPCON SCANDINAVIA A.B.

Neongatan 2; P.O.Box 25; 43151 Mölndal, SWEDEN
Phone: 46-(0)31-7109200 Fax: 46-(0)31-7109249 E-mail:medical@topcon.se; www.topcon.se

## TOPCON ESPAÑA S.A.

PRODE:+34-93-4734057 Fax:+34-93-4733932 E-mail: med PORTUGAL OFFICE Rua da Forte,6-6A,L-0.22,2790-072 Carnaxide,PORTUGAL Phone:+351-210-994626 Fax:+351-210-938786 www.topco

# TOPCON ( GREAT BRITAIN ) LTD.

use, Kennet Side,Bone Lane, Newbury,Berkshire RG14 5PX United Kingdom 0)1635-551120 Fax:+44-(0)1635-551170 E-mail:medical@topcon.co.uk; www.topcon.co.uk

# TOPCON POLSKA Sp. z. o. o. ul. Warszawska 23; 42-470 Siewierz, POLAND Phone: +48-(0)32-6705045 Fax: +48-(0)32-6715

# TOPCON SINGAPORE MEDICAL PTE. LTD.

lexandra Terrace, #08-27 The Comtech, Singapore 118502 ne:+65-68720606 Fax:+65-67736150 E-mail:medical\_sales@topcon.comm.sg www.topcon.com.sg

TOPCON INSTRUMENTS ( MALAYSIA ) SDN.BHD.
No. D1, (Ground Floor), Jalan Excella 2, Off Jalan Ampang Putra,
Taman Ampang Hilli, 55100 Kuala Lumpur.MALAYSIA
Phone:+80-(0)-342709866 Fax;+80-(0)-342709766

# TOPCON INSTRUMENTS (THAILAND) CO.,LTD.

77/162 Sinnsathorn Tower, 37th Floor, Krungthonburi Klongsarn, Bangkok 10600, THAILAND Phone:+66(0)2-440-1152~7 Fax:+66-(0)2-440-1158

# TOPCON CORPORATION BEIJING OFFICE

# TOPCON CORPORATION SHANGHAI OFFICE

14L Huamin Empire Plaza, No.726, Yan-an Xi Road, Shanghai, 200050, CHINA Phone:+86-(0)21-5238-7722 Fax:+86-(0)21-5237-0761

# TOPCON CORPORATION BEIRUT OFFICE

TOPCON CORPORATION DUBAI OFFICE P.O.Box 293705,Office C-25(row C-2), Dubai Airport Free Zone, Dubai, United Arad Emirates Phone:+971-4-299-5901

> PRINTED WITH SOY INK © 2010 TOPCON CORPORATION All Rights Reserved
> Printed in Japan 2010 06-20 NP 903-9 TMS • TS • AS